

● **About Us** ●

Here at Dr Physio USA, we make the most imaginative, unwinding, and restorative, well-being and individual care Products.

in the hubbub of today's lifestyle, it becomes impossible to nurture yourself at times. In the chaos of maintaining every aspect of life efficiently, we forget that our bodies are not machines, they need absolute care and tenderness. Pain in neck, shoulder, back or foot should never be taken for granted. They can lead to serious ailments in the long run. To help you ease the situations of painfulness, Dr Physio USA brings to you excellent ways of relaxation and pain relief through our premium gadgets. We offer a non-intrusive answer to help renew and reestablish your muscles and relieve throbbing pain with our various products

The variety of items and massagers we offer are intended to target particular areas of the body viz neck, shoulder, back, foot and so on. The adaptable capabilities of our gadgets are meant for utmost body relaxation after a long day at work.

Sit Back, Relax and enjoy the amazing Dr Physio experience by Dr Trust!!

Dr Physio from Dr Trust
www.drphys.io | www.drtrust.in

Dr Physio
Supervolt Go Massager - 1023 Pro


● **Quick Start Guide** ●

- | | |
|--|--|
| Step 1
Press the power button to start the message gun. | Step 2
Press the switch 1 to 6 times to change gears. |
| Step 3
Point the head at the muscle and move it along the muscle slowly. | Step 4
Long press the power button to turn the massager OFF. |

● **Overview** ●

Dr Physio Supervolt Go 1023 Pro is a handheld percussion massager that applies pulses of concentrated pressure deep into your muscle tissue. Deep-tissue massage works effectively on knots and tension through improving blood flow. Proper regular use of it can reduce your recovery time while relieving the pain, tension, and soreness that comes naturally from intense exercise.

● **Product Description** ●

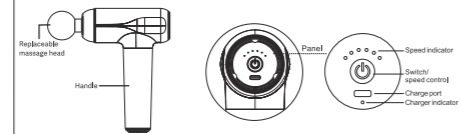

● **Massage Heads** ●

- Equipped with four replaceable message heads
- Round Massage Head**
The larger head is predominantly used for large muscle groups like the glutes, hamstrings, and quads, as it has a larger surface area.
- Bullet Massage Head**
This head is ideal for trigger point massage. It concentrates the pulses on a small surface area, so ideal for deep massage or to use on the hands and feet.


- Flat Massage Head**
The Flat Head has a firm, smooth, flat surface. It is very good for pressing into denser muscle groups, such as quads, glutes, calves, pecs. Again, it is a great multipurpose attachment.
- Fork Head**
With the fork, you can position the Achilles tendon in-between it, and the attachment will not put pressure on the tendon itself but treat the surrounding tissue.

● **Method Of Massage Head Installation** ●

Please turn off message gun and slowly insert the message head into the groove of the muzzle of it. When they are combined, the installation is finished.



Note:
It is forbidden to put your fingers, hair or other parts near the axis or back of the message head to avoid being pinched.

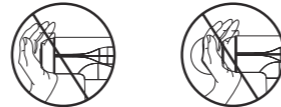

● **Operation Instructions** ●

- Step 1:**
Choose a message head you need and install it.
- Step 2:**
Press the power button to turn the message gun ON. By default, it enters into the first intensity level.

- Step 3:**
Press the button to increase the intensity of the massage.
- Step 4:**
It has 6 intensity levels and when it reaches the 6th level, press the button again to shut it off.
- Step 5:**
The massager will automatically shut down after 15 minutes.
- Step 6:**
Long press the power button to turn the massager OFF.

● **Safety Cautions** ●

Read below shared information before use, to avoid any damage or become dangerous!

- For use only by healthy adults other than older adults;
- Use the product only on the soft tissues of the body without causing pain or discomfort;
- Do not use on the head, or any hard or bony part of the body; Either controlled downward pressure or the application of pressure may produce bruising;
- Use only on dry, clean skin surfaces on body parts, pressing and moving gently on the skin;
- Keep your fingers, hair or other body parts away from the shaft and back of the vibrating head, as they may get caught;
- Do not apply to bony areas where there is less muscle coverage, such as the shins, backs of the feet, back of the hands;
- Do not use the massager on or around sensitive areas such as eyes, teeth, breast implants, genital organs, artificial limbs, etc;
- This product is not a medical device;
- Do not use it when you get bruised or feel painful and uncomfortable;
- Please check the main unit and cable before using, if damaged, it must be replaced by the manufacturer or its appointed professional to avoid danger;
- Check the massage area frequently and – stop using the product as soon as signs of pain or discomfort appear;

● **Warning** ●

- Avoid getting the message gun wet.
- Do not make any private modifications to the equipment.
- The use of this device on an injured part of the body is prohibited without consulting your doctor.

● **Disclaimer** ●

Not to be used by children or frail individuals. They should be accompanied by an adult when using any vibration device. Not to be used during pregnancy (definitely not on the abdomen area and you should take advice from your doctor if it is safe to use elsewhere on the body). Not to be used by cancer patients that are receiving treatment. Avoid massaging areas with localised skin infections, deep bruising or varicose veins. Not to be used if you have haemophilia, complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis or spondylolysis, recent joint replacements or IUDS, metal pins or plates or any concerns about your physical health.

Note: Not Intended for Therapeutic Use

These contraindications do not mean that you are not able to use this device but, we advise you to consult a doctor first.

● **Technical Specifications** ●

- Rated Voltage: 5V
- Rated Power: 25W
- Auto Off: 15 Minutes
- Weight: 0.61kg

● **Cleaning And Storage** ●

Wipe the surface of the device with a slightly damp towel and dry with a soft cloth.

● **Troubleshooting** ●

Please refer to the following tips and solutions if you encounter problems during use.

Problem	Possible Solution
Cannot work	No power.
Suddenly stop working	The message gun maybe about to be dead. Please use it after fully charged. Time is up. It will be auto off after 15 mins.
A slight noise when working	This is due to the sound generated by the mechanical structure, not a malfunction.

● **Box Content** ●

- 1 x Message Gun
- 1 x Cable
- 4 x Message Heads

● **Customer Support** ●

CONTACT ADDRESS
USA
Nureca INC.USA
276 5th Avenue, Suite 704-397,
New York (NY) - 10001, USA
INDIA
Corporate Office (Mumbai)
Nureca Limited
Office No. 101, 1st Floor, Udyog Bhavan
Sonawala Lane Goregaon East, Mumbai
Maharashtra - 400063, India
Contact us
India: +91-7527013265 / +91-9356658436
Website: www.drtrust.in
Corp Website: www.nureca.com
Email: customercare@nureca.com
Connect with us on social networks
Facebook: @drtrust
Instagram: @drtrustsin
Youtube: NurecaUsa
COPYRIGHT © 2024 NURECA LTD ALL RIGHTS RESERVED



Scan to View
Product Demo Video
www.drtrustusa.com/1023pro